

A.I.T. France's head office is located at Houilles near Paris. It occupies more than 1000 m² of workshop and offices, production and distribution facilities





DISTRIBUTED BY:

Greyhound Chromatography & Allied Chemicals

6 Kelvin Park, Birkenhead, Merseyside, CH41 1LT, UK

Tel: (+44) 0 151 649 4000 • Fax: (+44) 0 151 649 4001

Email: info@greyhoundchrom.com • Web: www.greyhoundchrom.com

Orders may be placed by phone, fax, email or web shop.

Since 1991, A·I·T· France has been an acknowledge chromatography specialist[.] Our legal and financial independence and our long experience have made us develop our wide range of smart caps :

The Smart Caps' range for more than 10 years, thus ranking us as one of the European leaders in this field, thanks to a workshop with 1000 m² of workshop[.] Our team is inspired by the will to offer to its clients a flexible and an efficient service with competitive priced products[.]

A·I·T· France is very attentive to your needs· Being an authentic manufacturer, we are ready to help you with particular requests and equipment installation in your laboratories, even for one single item without any further cost· A·I·T· France will always strive to earn your confidence·

Should you want to distribute our products ? Please, contact us, OEM manufacturing are welcome

Pascal AZNAR President of A·I·T· France SUMMARY

4-7

13 14

19

HEALTHY AND WASTE CAPS,

HOW IT WORKS ?

SMART HEALTHY CAPS

ANALYTICAL SYSTEM

Guide for installing a SHC	8
Smart Healthy Caps GL45	9
GL45 thread with stopcock	9-10
Ground neck	10
Distri Smart	11
Smart Healthy Caps with septum	12

PREPARATIVE HPLC SYSTEM

GL45 thread	
DIN 51/S55 & S60 thread	
Safety funnel for preparative	

ADAPTERS

Multi-Ways adapters	
Operating principle	
Classic adapters	

BOTTLE

DURAN[®] glass bottle GL45

ACCESSORIES

Air check valve 2 in 1 fittings Fittings and Ferrules Tubing



ACTIVATED CHARCOAL CARTRIDGE

Cartridge for Cans
Cartridge for Metallic barrel
Adapters
Adapters for JUSTRITE® can

SMART WASTE CAPS THREADS

_		
(Guide for installing a SWC	28
D	DIN 51 / S55	29
(GL45	30
В	Burkle S50	31
D	DIN 60 / S61	32
D	DIN 90 / S90	33
\land	VALGENE® bottle	34
	_evel indicator	35
\leq	Special demands	36

ADAPTERS

For cans	fitted	37-38

ACCESSORIES

15	2 in 1 universal fittings	39
	Leak tube fittings	40-41
	Conectors and adapters	42
16	Tubing and blind plugs	43
17	Quick connect fittings	44-45
18	- 3	

SAFETY FUNNELS

	Funnel with Hinged lid	48-49
20	Funnel with stopcock	50
	Funnel with ball valve	39-53
22	Electro conductive funnel	54
23		



56-57


Cans	56
Retention tray	57
Space-saving cans	57

NDEX

58-61

HEALTHY AND WASTE CAPS,

A ONE WAY AIR CHECK VALVE TECHNOLOGY

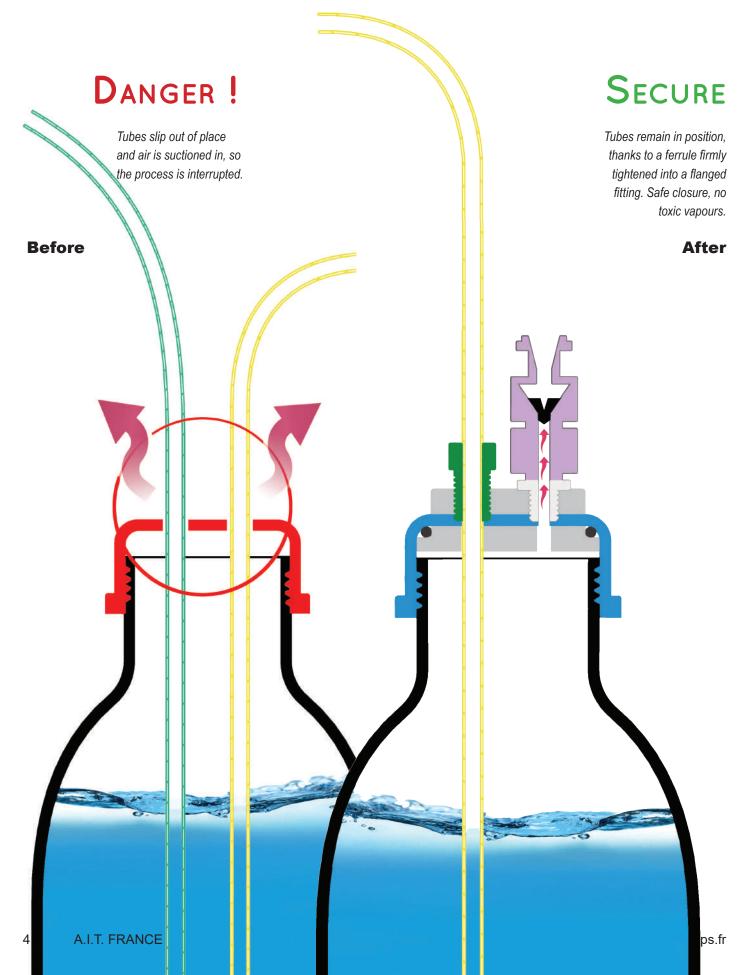

HOW DOES IT WORKS ?

A TECHNOLOGY TO PROTECT YOUR HEALTH 1/8" or 1/16" capillary - Air check valve Flanged fittings and ferrules -/ Body caps in Teflon (PTFE) Screw caps in polypropylene VITON[®] seal (PP) FRANCE http://www.smartcaps 5

HEALTHY AND WASTE CAPS,

A ONE WAY AIR CHECK VALVE TECHNOLOGY

Smart Healthy Caps are in response to modern laboratory safety regulations.

They are made of inert materials resistant to chemicals.


The use of a one-way air check valve allows for low cost safe and secure solvent extraction in your laboratory.

- ► No solvent evaporation. ✓
- Keep your concentration accurate for several days. \checkmark
- ► No mobile phase external contamination. ✓




The rotation of the central part of the caps allows a good ease of installation, and also avoids any entanglement of tubes.



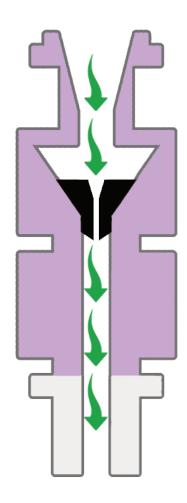

The membrane allows air to enter for pressure equalitization.

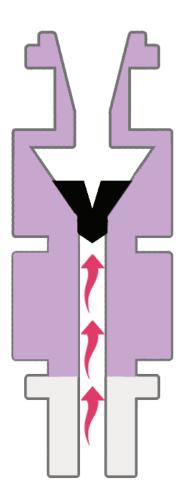


The membrane is a barrier against eluent vapours and gas leaks.

HOW DOES IT WORKS ?

A TECHNOLOGY TO PROTECT YOUR HEALTH


A Smart Waste Caps is formed by an airtight cap, drilled by a variable number of holes to host effluents tubes outgoing of the HPLC system.

The main piece is a cartridge filled with activated charcoal attached to the cap. This filter allows the absorption of 98% of the solvent's vapours released by the solvent's can.

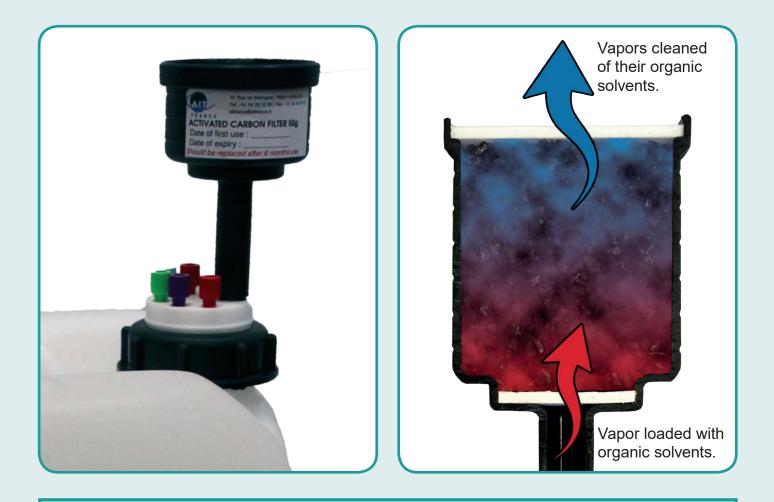


AIT France offers a diversified range of activated charcoal cartridge for the emission of eluent or solvent vapors; by special charcoal impregnation, the cartridges become an excellent barrier against concentrated acids or bases.

GUIDE FOR INSTALLING A SMART HEALTHY CAPS

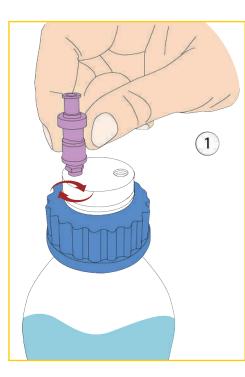


 Our Smart Healthy Caps are delivered with their preinstalled air check valves.
 However, if you need a valve that is resistant to polar solvent or amine derivatives vapors, simply unscrew the valve in place and then screw the new valve into its proper thread of the cap.

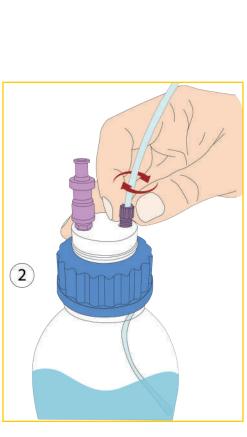

• Screw the 2 in 1 fitting into its proper thread on your Smart Healthy Cap. (See our "Installing a 2-in-1 fitting" guide on page 21 for more details)



• Installation completed!

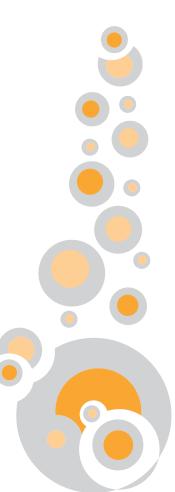



...And do not forget to indicate the date of commissioning of your air check valve in order to be able to anticipate its renewal!

ANALYTICAL SYSTEM



GL45 thread type with air check value

Smart Healthy Caps for GL45 or DURAN[®] types are the most widespread in laboratories. Smart Healthy Caps are delivered with their Universal 2 in 1 capillary fitting (external diameter 1/8" to 1/16").









DESCRIPTION	Ferrule *	Reference
Smart healthy caps GL45, 1 outlet and 1 air check valve	1/8"	SHC-1
Smart healthy caps GL45, 2 outlets and 1 air check valve	1/8"	SHC-2
Smart healthy caps GL45, 3 outlets and 1 air check valve	1/8"	SHC-3
Smart healthy caps GL45, 4 outlets and 1 air check valve	1/8"	SHC-4
Smart healthy caps GL45, 5 outlets and 1 air check valve	1/8"	SHC-5



* For another diameter, precise it when you order it.

Warning : Change your air check valve once a year.

GL45 THREAD WITH STOPCOCK.

GL45 Smart Healthy Caps with a stopcock offer more flexibility, without air bubbles inside the tubing. Smart Healthy Caps are delivered with their fittings 1/8" (external diameter 3,2 mm).



SHC-V2	DESCRIPTION	Ferrule *	Reference
	Smart healthy caps GL45, 1 outlet with shut-off and 1 air check valve	1/8"	SHC-V1
	Smart healthy caps GL45, 2 outlets with shut-off and 1 air check valve	1/8"	SHC-V2
shc-v4	Smart healthy caps GL45, 3 outlets with shut-off and 1 air check valve	1/8"	SHC-V3
	Smart healthy caps GL45, 4 outlets with shut-off and 1 air check valve	1/8"	SHC-V4
	*For another diameter, specify it when ordering.		
GROUND NECK	WITH AIR CHECK VALVE		7
	vith ground neck 29/32 are made for ground- Healthy Caps are delivered with their		

Universal 2 in 1 capillary fitting (external diameter 1/8" to 1/16").

SHC-1-EM



AD-45-EM

DESCRIPTION	FERRULE *	Reference
Smart healthy caps, ground neck (29/32 mm)	1 outlet (1/8"), 1 air check valve	SHC-1-EM
Smart healthy caps, ground neck (29/32 mm)	2 outlets (1/8"), 1 air check valve	SHC-2-EM
Smart healthy caps, ground neck (29/32 mm)	3 outlets (1/8"), 1 air check valve	SHC-3-EM
Smart healthy caps, ground neck (29/32 mm)	4 outlets (1/8"), 1 air check valve	SHC-4-EM

ANALYTICAL SYSTEM

DISTRI SMART SAFE WITH A SMART HEALTHY CAP

DESCRIPTION	Reference
Distri Smart Safe Adapter of 1 outlet (1/8" or 1/16") + 1 air check valve	SHC-1PV
Distri Smart Safe Adapter of 2 outlets (1/8" or 1/16") + 1 air check valve	SHC-2PV
Distri Smart Safe Adapter of 3 outlets (1/8" or 1/16") + 1 air check valve	SHC-3PV
Distri Smart Safe Adapter of 4 outlets (1/8" or 1/16") + 1 air check valve	SHC-4PV
Pipette 1 to 10 mL *	DIS 10
Pipette 2,5 to 25 mL *	DIS 25



ANALYTICAL SYSTEM

Smart Healthy Caps Septum offers you the possibility to introduce or withdraw a volume of solvent or reagent of your choice at any time.

DESCRIPTION	Reference
Smart Healthy Caps GL45 1 septum*, 1 outlet (1/8" or 1/16") and 1 air check valve	SHC-1-SEPT
Smart Healthy Caps GL45 1 septum*, 2 outlets (1/8" or 1/16") and 1 air check valve	SHC-2-SEPT
Smart Healthy Caps GL45 1 septum*, 3 outlets (1/8" or 1/16") and 1 air check valve	SHC-3-SEPT

*The septum is delivered as standard, non-pierced

SHC-2-SEPT

PREPARATIVE SYSTEM

SMART HEALTHY CAPS FOR GL45 SCHOTT® OR DURAN® TYPE



DESCRIPTION	REFERENCE
Smart Healthy Caps GL45 for preparative HPLC with 1 outlet (1/4") + 1 air check valve	SHC-P-A
Smart Healthy Caps GL45 for preparative HPLC with 2 outlets (1/4" and 1/8") + 1 air check valve	SHC-P-A-C
Smart Healthy Caps GL45 for preparative HPLC with 1 outlet (3/16") + 1 air check valve	SHC-P-B
Smart Healthy Caps GL45 for preparative HPLC with 2 outlets (3/16") + 1 air check valve	SHC-P-2B



SHC-P-A-C



SHC-P-A



SHC-P-B

OTHER CONFIGURATIONS POSSIBLE ON REQUEST!

PREPARATIVE SYSTEM

SMART HEALTHY CAPS FOR DIN51/S55 AND DIN60/S60

THREADS


Specially adapted for large volumes of preparative HPLC eluents. The air check valve works with a flow rate from 150ml / min to 500ml / min.



DESCRIPTION*	Reference
Smart Healthy Caps S55 for preparative HPLC with 1 outlet (1/4") + 1 air check valve	SHC55-P-A
Smart Healthy Caps S55 for preparative HPLC with 2 outlets (1/4" et 1/8") + 1 air check valve	SHC55-P-A-C
Smart Healthy Caps S55 for preparative HPLC with 1 outlet (3/16") + 1 air check valve	SHC55-P-B
Smart Healthy Caps S60 for preparative HPLC with 1 outlet (1/4") + 1 air check valve	SHC60-P-A
Smart Healthy Caps S60 for preparative HPLC with 2 outlets (1/4" et 1/8") + 1 air check valve	SHC60-P-A-C
Smart Healthy Caps S60 for preparative HPLC with 1 outlet (3/16") + 1 air check valve	SHC60-P-B



*also available with any type of thread (Nalgene, DIN90...)

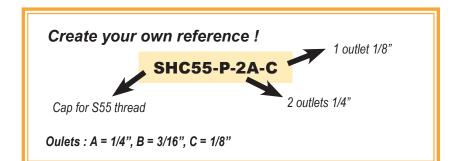

PREPARATIVE SYSTEM

SMART HEALTHY CAPS WITH FUNNEL FOR PREPARATIVE HPLC BOTTLE


This Smart Healthy Caps, equipped with a funnel with hinged lid, allows to pour the eluent into the bottle without having to unscrew the cap. It is equipped with an air check valve.

	DESCRIPTION	THREAD	Reference
		GL45	SHC45-FSP-B
	Funnel with hinged lid for	S55	SHC55-FSP-B
SHC45-FSP-B	solvent reservoir with one outlet of 4 mm (3/16") and	S60	SHC60-FSP-B
	one air check valve	NALGENE [®] B53	SHCB53-FSP-B
		Metallic barrel	SHCBA-FSP-B
		GL45	SHC45-FSP-A
le il	Funnel with hinged lid for	S55	SHC55-FSP-A
	solvent reservoir with one outlet of 6.4 mm (1/4") and	S60	SHC60-FSP-A
	one air check valve	NALGENE [®] B53	SHCB53-FSP-A
		Metallic barrel	SHCBA-FSP-A
		GL45	SHC45-FSP-C
	Funnel with hinged lid for	S55	SHC55-FSP-C
	solvent reservoir with one outlet of 3.2 mm (1/8") and	S60	SHC60-FSP-C
	one air check valve	NALGENE [®] B53	SHCB53-FSP-C
		Metallic barrel	SHCBA-FSP-C
	Thread type	<pre>c own reference ! SHC55 - FSP - B 4", B = 3/16", C = 1/8"</pre>	Tube outer diameter
http://www.smartcaps.fr			A.I.T. FRANCE 15

ADAPTERS WITH MULTIPLE WAYS

Multi-ways adapter allowing to connect a bottle with a GL45 thread to 2 & 4 GL45 safety caps*

The multi-ways adapter allows several channels to be purged in the eluent bottle while remaining airtight.

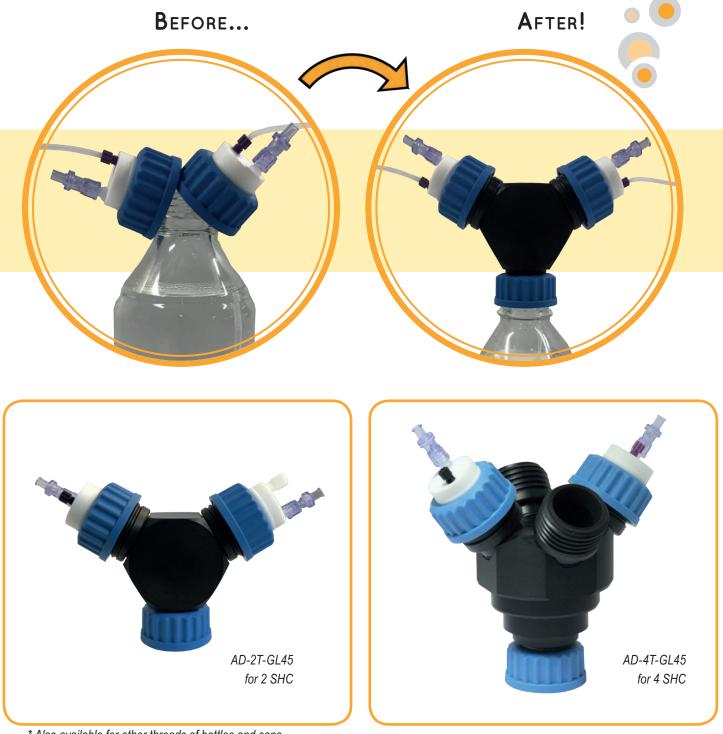

ADAPTERS

OPERATING PRINCIPLE

We can supply you with all types of adapters not described in the catalog. Do not hesitate to ask us!



32 mm

GL32 bottle Measures 32 mm for external threads. Smart Healthy Caps GL45



38 mm

GL38 bottle Measures 38 mm for external threads.







45 mm

GL45 bottle Measures 45 mm for external threads.

CLASSIC ADAPTERS

To connect a GL45 cap on a bottle with a different thread.

THREAD*	То**	MATERIAL	REFERENCE
GL40	GL45	PTFE	AD-40-45-T
GL40	GL45	PP	AD-40-45-P
GL32	GL45	PTFE	AD-32-45-T
GL32	GL45	PP	AD-32-45-P
GL38	GL45	PTFE	AD-38-45-T
GL38	GL45	PP	AD-38-45-P

AD-40-38-P

AD-40-45-P

PTFE : Teflon.

PP : Polypropylene.

* Bottle thread

* * Cap's Thread : Smart Healthy Caps.

Notes

•••••	• • • • • • • • • • • • • • • • • • • •	• • • • • • • • • • • • • • • • • • •	• • • • • • • • • • • • • • • • •
•••••		•••••	• • • • • • • • • • • • • • • • • • •
•••••		•••••	• • • • • • • • • • • • • • • • •
• • • • • • • • • • • • • •		•••••	• • • • • • • • • • • • • • • • •
• • • • • • • • • • • • • •		•••••	• • • • • • • • • • • • • • • • •
•••••		•••••	• • • • • • • • • • • • • • • • • •

BOTTLE

DURAN[®] GL45 BOTTLE

A complete range of solvent bottles perfectly adapted for GL45 Smart Healthy Caps.

VOLUME (ML)	GLASS	REFERENCE
250	Clear	GL250-C
250	Amber	GL250-A
500	Clear	GL500-C
500	Amber	GL500-A
1000	Clear	GL1000-C
1000	Amber	GL1000-A
2000	Clear	GL2000-C
2000	Amber	GL2000-A
5000	Clear	GL5000-C
5000	Amber	GL5000-A



400

The Code of

Accessories : Air check valve

All our Smart Healthy Caps are equipped with a standard air check valve.

For optimal use, it is recommended to *change the air check valve once a year!*

STANDARD AIR CHECK VALVE

For organo-aqueous eluents.

Weak and strong acidic eluents.

Non-amino basic eluents.

Polar & apolar eluents.

Maximum recommended flow: 100mL / min

REFERENCE	QTY
SHC-SV	1
SHC-SV10	10
SHC-SV50	50





AIR CHECK VALVE FOR PREPARATIVE HPLC

Same characteristics as the standard air check valve but for HPLC systems with very high flow rates (>100mL / min and <500mL / min).

AIR CHECK VALVE FOR BASIC AMINE SOLVENTS


Same characteristics as the standard air check valve, but essential for the use of eluents containing bases with high pKa (triethylamine) or high concentrations (ammonia).

Reference	QTY
SHC-SV-AB	1
SHC-SV-AB10	10
SHC-SV-AB50	50


Accessories : the 2 in 1 universal fitting

- Facilitates the connection of capillary tubes from 1 mm to 3,2 mm (1/8")
- The 2 in 1 Universal fitting replaces the entire fitting + ferrule 1/8" and 1/16"
- The 2 in 1 Universal fitting ensures excellent sealing


Guide for 2 in 1 fitting change:

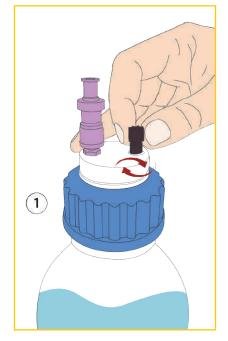

• Insert the 2 in 1 fitting into its proper thread, without fully tightening it.

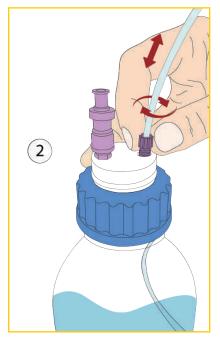

• Insert your tube into the 2 in 1 fitting, adjust its length and then screw the fitting until the tube stops moving.

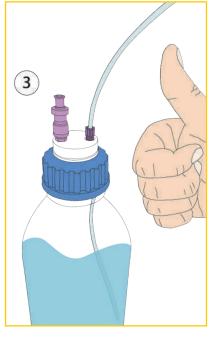

Installation completed!



DESCRIPTION	For Tube (OD)	Color	Q τγ/Ρκ	REFERENCE
2 in 1 Fitting, Type I Black	for 1/16" to 1mm	Black	5	CC-ONE-1BL
2 in 1 Fitting, Type I Purple	for 1/16" to 1mm	Purple	5	CC-ONE-1V
2 in 1 Fitting, Type II Black	for 1/8" to 1/16"	Black	5	CC-ONE-2BL
2 in 1 Fitting, Type II Purple	for 1/8" to 1/16"	Purple	5	CC-ONE-2V

Accessories

NUTS AND FERRULES FOR 1/8" AND 1/16"

Ferrules 1/8", ref : F18. Flangefree fitting for ferrule 1/8" OD. Ferrules 1/16", ref : F16. Flangefree fitting for ferrule 1/16" OD.



* Please join to the reference the color you want. For example : R18-N for black nuts.		Тиві Імсн	e (OD) Millimeter	Reference Nuts per 10 units *	REFERENCE FERRULES PER 10 UNITS
		1/4"	6,3 mm	R14	F14
nato.		3/16"	4 mm	R316	F316
Red : R Blue: B		1/8"	3,2 mm	R18	F18
Green: G			2,5 mm	R18	F25
Purple: P Black: N		1/16"	1,6 mm	R16	F16
			1,0 mm	R16	F10

BLIND PLUG, ALLOWING YOU TO CLOSE AN UNUSED OUTLET OF YOUR SMART HEALTHY CAPS



PTFE TUBES



OD	ID	REFERENCE *	QUANTITY
1/4"	6,3 mm	TP14	1 mètre
3/16"	3 mm	TP316	1 mètre
1/8"	1,6 mm	TP18	1 mètre
1/16"	1,0 mm	TP16	1 mètre



OD	ID	Reference *	QUANTITY
9,5 mm	6,4 mm	TYG-64	1 mètre
19,1 mm	12,7 mm	TYG-127	1 mètre



Notes

• • • • • • • • • • • • • • • • • • • •	• • • • • • • • • • • • • • • • • • • •	• • • • • • • • • • • • • • • • • • • •
•••••	•••••	• • • • • • • • • • • • • • • • • • • •
•••••	•••••	• • • • • • • • • • • • • • • • • • • •
• • • • • • • • • • • • • • • • • • • •	• • • • • • • • • • • • • • • • • • • •	
•••••	•••••	
•••••	•••••	• • • • • • • • • • • • • • • • • • • •
• • • • • • • • • • • • • • • • • • • •	•••••	• • • • • • • • • • • • • • • • • • • •

SMART WASTE CAPS ACTIVATED CHARCOAL CARTRIDGE ACTIVATED CHARCOAL CARTRIDGE FOR CANS FROM 2,5 TO 20 LITERS VATED CARBON FILTE N FILTER 100 VATED CARBON FILTER 509 Tel 01 34 93 10 85 Fm , 51 54 AIT FILTRE A CHARBON ACTIF 100g FC-25 FC-50 FC-100 Date de mise en serv Date de péremption : Each cartridge is delivered with a label on which you can write the commissioning date and the replacement date of your cartridge. The replacement date depends on the quantity of activated charcoal in the cartridge. * For each charcoal filter cartridge order, a time monitoring is delivered if you ask for. Just add the letter T at the end of the reference.

DESCRIPTION	QUANTITY OF CHARCOAL	DURATION	Reference
Activated charcoal cartridge	25 gr	3 months	FC-25
Activated charcoal cartridge	50 gr	6 months	FC-50
Activated charcoal cartridge with a lid	50 gr	6 months	FC-50-C
Activated charcoal cartridge	100 gr	12 months	FC-100
Activated charcoal cartridge loaded with H ₂ SO ₄ for volatile base in waste eluants.	100 gr*	6 months	FC-100-H2SO4
Activated charcoal cartridge loaded with KC for volatile acid in waste eluants.	0H 100 gr*	6 months	FC-100-KOH

Cartridge of activated charcoal for metallic barrel from 50 to 200 liters

For a metallic barrel, a higher quantity of activated charcoal is needed. We propose you 2 models of filter depending on the volume of your metallic barrel.

DESCRIPTION	FOR A METALLIC BARREL OF	THREAD	REFERENCE
Activated charcoal cartridge of 250 gr	50 L	2" Mausser	FC-250
Activated charcoal cartridge of 500 gr	200 L	2" Mausser	FC-500


A wide range of safety funnels is proposed to you page 46 to 55.



ACTIVATED CHARCOAL CARTRIDGE

Adapters for activated charcoal cartridge

If you have a problem to place the can with the cartridge under the bench, we offer you 4 types of adapters which can be combined to solve all of you problems. You can see below few types of combinations. Other combinations are possible.

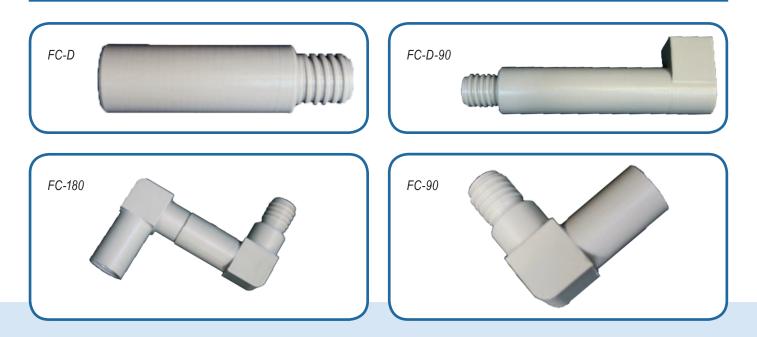

DESCRIPTION	Reference	
Straight extension long piece	FC-D	FC-90
90° Adapter for angled connection	FC-90	
Double angled connection adapter (combination of 2 FC-90)	FC-180	
Straight extension long piece with 90° adapter for angled connection	FC-D-90	



http://www.smartcaps.fr

SPECIAL ADAPTERS* FOR JUSTRITE® CAN


* Allows to accept activated charcoal cartridge on JUSTRITE[®] can.




Adapter with 3 entries and a tube fitting for a direct CPC connection	AD-CPC
Adapter with 3 entries and a tube fitting for a connection on $JUSTRITE^{ extsf{s}}$ can	AD-D-JUST
Adapter with an activated charcoal cartridge location for JUSTRITE [®] container	FC-AD-JUST
Adapter with an activated charcoal cartridge location for CPC connection	FC-AD-CPC

Smart Waste Caps

Guide for installing a Smart Waste Caps

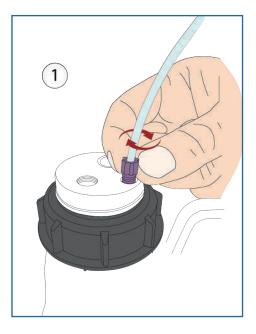

Insert the 2 in 1 fitting into its proper thread

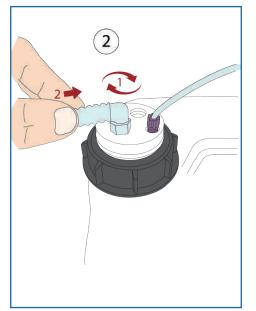


 Screw the Leak fitting into its proper thread on your Smart Waste Cap.
 Connect your tube to the Leak fitting


HEY!

DON'T

FORGET!

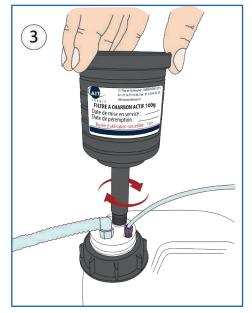

• Screw your activated charcoal cartridge into its proper thread on your Smart Waste Caps

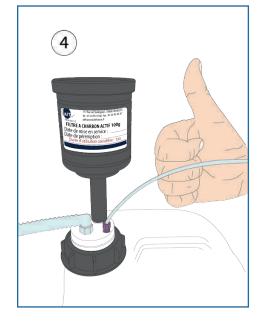

Installation completed!

And do not forget to indicate the date of commissioning and the expiration date of your activated charcoal cartridge to anticipate its renewal!

SMART WASTE CAPS FOR DIN 51 / S55

AIT Smart Waste Caps for DIN 51/S55 with their Universal 2 in 1 capillary fitting (external diameter 1/8" to 1/16") and tube fitting 6 to 9 mm.



SW55-4-C : 4 entries 1/8" and charcoal cartridge emplacement




SW55-2-C-1L : 2 entries 1/8", 2 tube fittings and charcoal cartridge emplacement



SW55-1-C-3L : 1 entry 1/8", 3 tube fittings and charcoal cartridge emplacement

If you can't find the reference you want, then, create your own Smart Waste Cap even just for one item, without further cost.

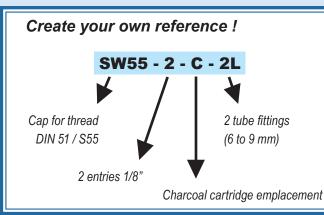





SW55-4-C-1L : 4 entries 1/8", a charcoal cartridge emplacement and 1 tube fitting



SW55-5-C : 5 entries 1/8" and a charcoal cartridge emplacement



SW55-3-C-2L : 3 entries 1/8", a charcoal cartridge emplacement and 2 tube fittings



SW55-3-C-3L : 3 entries 1/8", a charcoal cartridge emplacement and 3 tube fittings

SMART WASTE CAPS THREADS

SMART WASTE CAPS FOR GL45

AIT Smart Waste Caps for GL45 with their Universal 2 in 1 capillary fitting (outer diameter 1/8" to 1/16") and tube fitting.

SW45-2-C : 2 entries 1/8" and charcoal cartridge emplacement



SW45-5-C : 5 entries 1/8" and charcoal cartridge emplacement




SW45-4-C-1L : 4 entries 1/8", 1 tube fitting and charcoal cartridge emplacement



SW45-1-C-2L : 1 entry 1/8", 2 tube fittings and charcoal cartridge emplacement

If you can't find the reference you want, then, create your own Smart Waste cap even just for one item, without further cost.

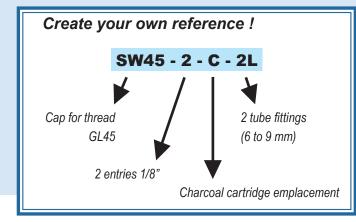







SW45-C-2L : a charcoal cartridge emplacement and 2 tube fittings 30 A.I.T. FRANCE

SW45-1-C-1L : 1 entry 1/8", a charcoal cartridge emplacement and 1 tube fitting



SW45-2-C-2L : 2 entries 1/8", a charcoalcartridge emplacement and 2 tube fittings



SW45-C-3L : a charcoal cartridge emplacement and 3 tube fittings

http://www.smartcaps.fr

Smart Waste Caps threads

SMART WASTE CAPS BURKLE AND DIN50

Smart Waste Caps for space-saving cans S50 equipped with their Universal 2 in 1 capillary fitting (outer diameter 1/8" to 1/16") and tube fitting 6 to 9 mm.



SWB-1-C-1L : 1 entry 1/8", 1 tube fitting and charcoal cartridge emplacement



SWB-3-C : 3 entries 1/8" and charcoal cartridge emplacement




SWB-1-C-3L : 1 entry 1/8", 3 tube fittings and charcoal cartridge emplacement



SWB-5-C-1L : 5 entries 1/8", 1 tube fitting and charcoal cartridge emplacement

If you can't find the reference you want, then, create your own Smart Waste cap even just for one item, without further cost.

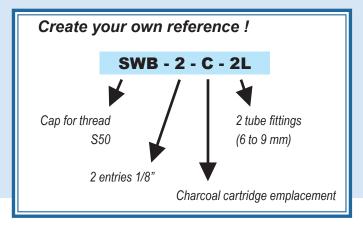



SWB-5-C: 5 entries 1/8" and charcoal cartridge emplacement

SWB-2-C-1L : 2 entries 1/8", 1 tube fitting and charcoal cartridge emplacement



SMART WASTE CAPS THREADS

SMART WASTE CAPS FOR DIN 60 / S61

AIT Smart Waste Caps for DIN 60/S61 with their Universal 2 in 1 capillary fitting (outer diameter 1/8" to 1/16") and tube fitting 6 to 9 mm.

SW60-2-C : 2 entries 1/8" and charcoal cartridge emplacement



SW60-2-C-1L : 2 entries 1/8", 1 tube fitting and charcoal cartridge emplacement




SW60-5-C-1L : 5 entries 1/8", 1 tube fitting and charcoal cartridge emplacement



SW60-1-C-4L : 1 entry 1/8", 4 tube fittings and charcoal cartridge emplacement

If you can't find the reference you want, then, create your own Smart Waste cap even just for one item, without further cost.

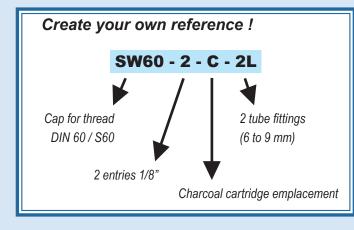





SW60-1-C-1L : 1 entry 1/8", a charcoal cartridge emplacement and 1 tube fitting 32 A.I.T. FRANCE



SW60-6-C : 6 entries 1/8" and a charcoal cartridge emplacement



SW60-2-C-4L : 2 entries 1/8", a charcoal cartridge emplacement and 4 tube fittings



SW60-C-1L : a charcoal cartridge emplacement and 1 tube fitting

http://www.smartcaps.fr

SMART WASTE CAPS THREADS

SMART WASTE CAPS FOR DIN 90 / S90

AIT Smart Waste Caps for DIN 90/S90 with their Universal 2 in 1 capillary fitting (outer diameter 1/8" to 1/16") and tube fitting 6 to 9 mm.

SW90-3-C : 3 entries 1/8" and charcoal cartridge emplacement



SW90-3-C-1L : 3 entries 1/8", 1 tube fitting and charcoal cartridge emplacement




SW90-3-C-2L : 3 entries 1/8", 2 tube fittings and charcoal cartridge emplacement



SW90-4-C-1L : 4 entries 1/8", 1 tube fitting and charcoal cartridge emplacement

If you can't find the reference you want, then, create your own Smart Waste cap even just for one item, without further cost.

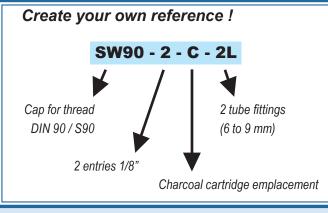





SW90-4-C : 4 entries 1/8" and a charcoal cartridge emplacement



SW90-2-C : 2 entries 1/8" and a charcoal cartridge emplacement



SW90-2-C-1L : 2 entries 1/8", a charcoal cartridge emplacement and 1 tube fitting



SW90-2-C-2L : 2 entries 1/8", a charcoal cartridge emplacement and 2 tube fittings A.I.T. FRANCE 33

SMART WASTE CAPS THREADS

SMART WASTE CAPS FOR NALGENE'S® BOTTLES

AIT Smart Waste Caps for NALGENE's[®] bottles with their Universal 2 in 1 capillary fitting (outer diameter 1/8" to 1/16") and tube fitting 6 to 9 mm.

	THREAD	Reference		THREAD	Reference
2	B53	SWB53-C-1L		B53	SWB53-C-3L
	B83	SWB83-C-1L		B83	SWB83-C-3L
SWB83-C-1L			SWB83-C-3L		
	THREAD	Reference		THREAD	Reference
	Thread B53	REFERENCE SWB53-2-C			REFERENCE SWB53-4-C-1L
SWB83-2-C			SWB83-4-C-1L	B53	

For the tread of the tread

SMART WASTE CAPS THREADS

LEVEL INDICATOR

Thanks to the level indicator and its integrated floater, you can check immediately the can level.

The level indicator is absolutely necessary for all the opaque cans and the metallic barrel. We can also equipped the level indicator with fittings 1/8" or 1/16" and tube fitting 6 to 9 mm.

How it works ?

The level indicator is an option. So you can add it to the Smart Waste Caps reference that you want to equip.

For example, you want to equip the Smart Waste Cap following : SW45-2-C-1L with a level indicator. It becomes :

SW45-2-C-1L-LI

Please do not hesitate to contact us with any questions you may have.



SPECIAL DEMANDS

These 2 elements are easily fixed via their screw system and can go through the bench or be at the edge of the bench.



Through-hole adapter for SW45 caps with 1 outlet for tube fitting (6 to 9mm in standard)

Multi-inlet system,6 1/8" entries and 4 entries for tubes, with an outlet to a can under the bench (11 to 13mm in standard)

Special manufacturing at the request of any cap or accessory according to your own specifications: from a single piece to more than a hundred!

ADAPTERS

ALL TYPES OF ADAPTERS ARE POSSIBLE. UNLIMITED COMBINATIONS! (AS FAR AS POSSIBLE)



45 mm

GL45 bottle Measures 45 mm for external thread.



48 mm

S55 can

Measures 48 mm for external thread.

SW45-C-1L -



Adapter 53 mm S60 can Measures 53 mm for external thread. Thread -

FILTRE A CHARBON ACTIF 100g

d'utilisation conseill


Date de mise en service : Date de péremption :

ADAPTERS

PP OR PTFE ADAPTERS FOR CANS EQUIPPED WITH A SMART WASTE CAPS



{	Male
)	Female

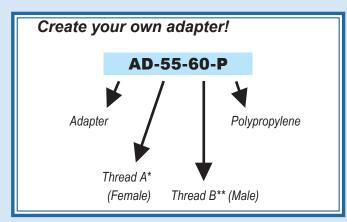

THREAD A* (female)	THREAD B** (MALE)	Reference ***
GL45	GL32	AD-45-32-P
GL40	GL32	AD-40-32-P
GL45	GL38	AD-45-38-P
GL40	GL38	AD-40-38-P
GL40	GL45	AD-40-45-P
GL45	S55	AD-45-55-P
S55	S60	AD-55-60-P
S55	S90	AD-55-90-P
S60	S90	AD-60-90-P

***also available in Teflon: for this, replace "P" with "T" in the reference of the desired product


*A = can thread

**B = Smart Waste Caps thread


Accessories : 2 in 1 Universal fitting

- Facilitates the connection of capillary tubes from 1 mm to 3,2 mm (1/8")
- The 2 in 1 Universal fitting replaces the entire nut + ferrule 1/8" et 1/16"
- The 2 in 1 Universal fitting ensures excellent sealing


Guide to change a 2 in 1 fitting :

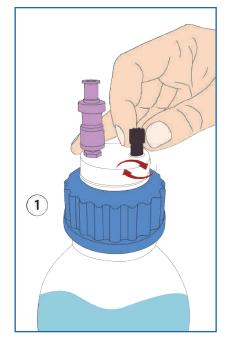

• Insert the 2 in 1 fitting into its proper thread, without fully tightening it.

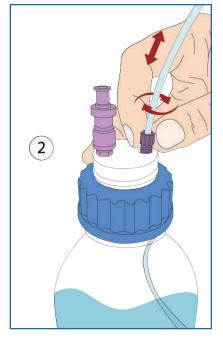


• Insert your tube into the 2 in 1 fitting, adjust its length and then screw the fitting until the tube stops moving.

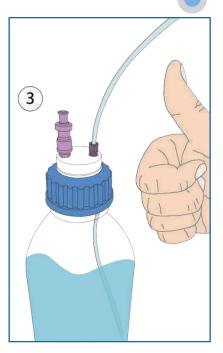

Installation completed!

DESCRIPTION	For tubes (OD)	Color	QTY/PK	REFERENCE
Fitting 2 in 1 Type I Black	of 1/16" to 1mm	Black	5	CC-ONE-1BL
Fitting 2 in 1 Type I Violet	of 1/16" to 1mm	Purple	5	CC-ONE-1V
Fitting 2 in 1 Type II Black	of 1/8" to 1/16"	Black	5	CC-ONE-2BL
Fitting 2 in 1 Type II Violet	of 1/8" to 1/16"	Purple	5	CC-ONE-2V



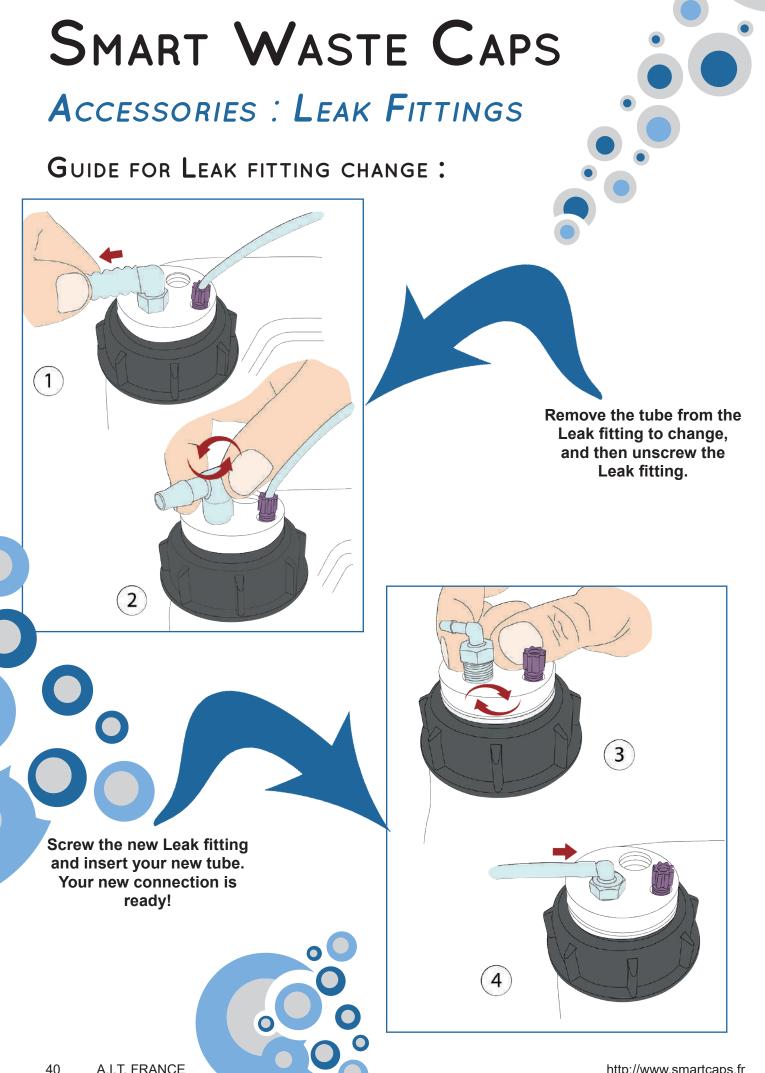


A.I.T. FRANCE

Accessories : leak fittings*

*The illustrations of the products are in CAD and are not contractual. The leak fittings are made of translucent PP.



STRAIGHT LEAK FITTING / PER UNIT

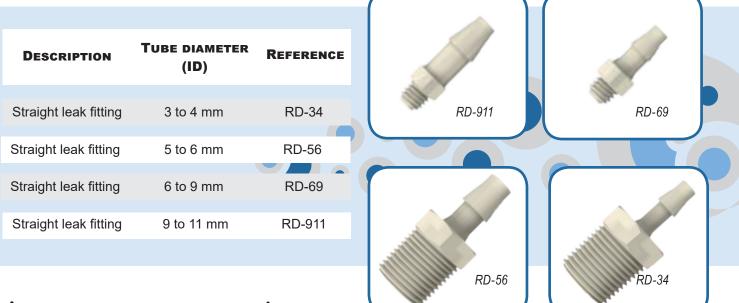

ANGLED LEAK FITTING / PER UNIT

DESCRIPTION	Tube diameter (ID)	Reference		
Angled leak fitting	3 to 4 mm	RC-34		
Angled leak fitting	5 to 6 mm	RC-56	RC-56	RC-69
Angled leak fitting	6 to 9 mm	RC-69		
Angled leak fitting	9 to 11mm	RC-911		

DOUBLE ANGLED LEAK FITTING / PER UNIT

DESCRIPTION	ID	REFERENCE
Double angled leak fitting	5 to 6 mm	RCY-56
Double angled leak fitting	6 to 9 mm	RCY-69




ACCESSORIES : CONNECTORS AND ADAPTERS



Y CONNECTORS*, PER UNIT

DESCRIPTION	TUBE ID	REFERENCE
Y-Connector 3 mm	3 mm	Y-3
Y-Connector 5 mm	5 mm	Y-5
Y-Connector 5 mm	5 mm	G-1
Y-Connector 7 mm	7 mm	Y-7
Y-Connector 9 mm	9 mm	Y-9
Y-Connector 11 mm	11 mm	Y-11

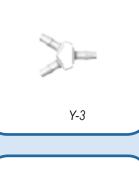








CONICAL REDUCERS AND ADAPTERS

For different sizes of tubes

INLET DIAMETER	OULET DIAMETER	Reference
4,2 - 8 mm	4,2 - 8 mm	CR-1
10,5 - 16 mm	6,5 - 12 mm	CR-3
11 - 15 mm	7 - 12 mm	CR-6

*The illustrations of the products are in CAD and are not contractual. The connectors are made of translucent PP.

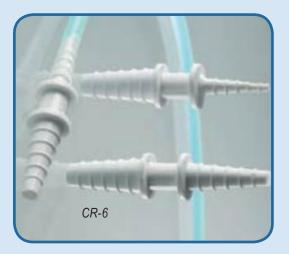

SMART WASTE CAPS Accessories : TUBES AND BLIND PLUG

PTFE TUBES



0	D	ID REF	ERENCE* G	QUANTITY
1/4	4" 6,3	3 mm	TP14	1 meter
3/1	6" 3	mm T	ГР316	1 meter
1/8	3" 1,6	6 mm	TP18	1 meter
1/1	6" 1,0) mm	TP16	1 meter

TYGONTM TUBES

OD	ID	Reference *	QUANTITY
9,5 mm	6,4 mm	TYG-64	1 meter
19,1 mm	12,7 mm	TYG-127	1 meter



BLIND PLUGS

THEY CLOSE THE UNUSED INLETS OF YOUR SMART WASTE CAPS

		\bigcirc		For	REFERENCE	QUANTITY
		-	T	1/8" or 1/16" fitting	VP-F18	Per 10
		11		Leak fitting	VP-LEAK	Per 5
Ĩ				Charcoal cartridge fitting	VP-CA	Unit
V	P-F18	VP-LEAK	VP-CA			

SMART WASTE CAPS Accessories : Quick Connect Fittings



Female Part

DESCRIPTION	OUTLET	PASSAGE	Reference	Α	В
	1/16" (1,6 mm)	3,2mm	S	22,35	35,56
Quick Connect (female	1/8" (3,2 mm)	3,2mm	Т	22,35	41,91
part) without shutoff valve	1/4" (6,4 mm)	3,2mm	U	22,35	46,99
	3/8" (9,5mm)	6,4mm	V	25,40	49,53
Quick Connect (female part) with shutoff valve	1/16" (1,6 mm)	3,2mm	W	22,35	35,56
	1/8" (3,2 mm)	3,2mm	Х	22,35	41,91
	1/4" (6,4 mm)	3,2mm	Y	22,35	46,99
	3/8" (9,5mm)	6,4mm	Z	25,40	49,53

ASK US FOR ADVICE FOR YOUR CUSTOM CONFIGURATIONS!

SMART WASTE CAPS Accessories : Quick Connect Fittings

CHOOSE THE MALE QUICK CONNECT FITTING THAT BEST SUITS YOUR NEEDS AND MODULATE THE CONNECTION ACCORDING TO YOUR WISHES. SW55-C-Z & QC64-38-CD

MALE PART

-
-
-
22,86
-

-	In-Line	1/4" (6,4mm)	6,4mm	QC64-14-CD-C	QC64-14-CD	16,00	34,29/ 50,80	-
and the	Hose Barb	3/8" (9,5mm)	6,4mm	QC64-38-CD-C	QC64-38-CD	16,00	34,29/ 46,23	-
)	Elbow	1/4" (6,4mm)	6,4mm	QC64-14-CC-C	QC64-14-CC	16,00	30,99/ 34,54	24,38
J.	Hose Barb	3/8" (9,5mm)	6,4mm	QC64-38-CC-C	QC64-38-CC	16,00	32,51/ 36,06	24,38

ATTENTION: MALE / FEMALE COMPATIBILITY IS ONLY ENSURED IF THE PASSAGE DIAMETER IS THE SAME FOR BOTH PARTS!

A wide range to meet all your needs











FUNNEL WITH HINGED LID FOR RESIDUAL SOLVENTS

It is installed on a can or a barrel of residual solvent. It is practical and you can equip it with an activated charcoal cartridge if the storage is in a non-ventilated room.

The activated charcoal cartridge prevents the creation of an overpressure in the can if it is exposed to heat.



DESCRIPTION	REFERENCE WITH CARTRIDGE [*]	REFERENCE WIT- HOUT CARTRIDGE
Smart Waste Caps funnel with hinged lid for can S55	SW55-FSC	SW55-FSB
Smart Waste Caps funnel with hinged lid for can S60	SW60-FSC	SW60-FSB
Smart Waste Caps funnel with hinged lid for can S90	SW90-FSC	SW90-FSB
Smart waste caps funnel with hinged lid for space-saving can S	50 SWB-FSC	SWB-FSB
Smart Waste Caps funnel with hinged lid for metallic barrel	SWBA-FSC	SWBA-FSB


Funnel with hinged lid for residual solvents

This funnel has a double interest: the direct discharge through the funnel of all your waste, but also the collection of effluent from your HPLC systems with the integration of inlet fittings and / or leak fittings on the funnel.

Funnel with hinged lid with 1 activated charcoal cartridge, 2 inlets 1/8" (3.2mm) and 2 leak fittings (6 to 9mm) I Sp Funnel with hinged lid with 1 activated charcoal cartridge,	THREAD S55 S60 S90 Metallic barrel S55 S60 S90 Metallic barrel S90 S90 Metallic barrel S90 S90<	REFERENCE SW55-2-1L-FSC SW60-2-1L-FSC SW90-2-1L-FSC SWBA-2-1L-FSC SWB-2-1L-FSC SW60-2-2L-FSC SW90-2-2L-FSC SWBA-2-2L-FSC SWBA-2-2L-FSC SWBA-2-2L-FSC SWBA-2-2L-FSC	SW55-2-1L-FSC
Activated charcoal cartridge, 2 inlets 1/8" (3.2mm) and 1 leak fitting (6 to 9mm) Funnel with hinged lid with 1 activated charcoal cartridge, 2 inlets 1/8" (3.2mm) and 2 leak fittings (6 to 9mm) Funnel with hinged lid with 1 activated charcoal cartridge,	S60 S90 Metallic barrel S55 S55 S60 S90 Metallic barrel S90 Metallic barrel	SW60-2-1L-FSC SW90-2-1L-FSC SWBA-2-1L-FSC SWB-2-1L-FSC SW55-2-2L-FSC SW60-2-2L-FSC SW90-2-2L-FSC SWBA-2-2L-FSC SWBA-2-2L-FSC SWBA-2-2L-FSC SWBA-2-2L-FSC	SW55-2-1L-FSC
Activated charcoal cartridge, 2 inlets 1/8" (3.2mm) and 1 leak fitting (6 to 9mm) Funnel with hinged lid with 1 activated charcoal cartridge, 2 inlets 1/8" (3.2mm) and 2 leak fittings (6 to 9mm) Funnel with hinged lid with 1 activated charcoal cartridge,	S90 Metallic barrel pace-saving can S55 S60 S90 Metallic barrel pace-saving can	SW90-2-1L-FSC SWBA-2-1L-FSC SWB-2-1L-FSC SW55-2-2L-FSC SW60-2-2L-FSC SW90-2-2L-FSC SWBA-2-2L-FSC SWBA-2-2L-FSC	SW55-2-1L-FSC
2 inlets 1/8" (3.2mm) and 1 leak fitting (6 to 9mm)	Metallic barrel space-saving can S55 S60 S90 Metallic barrel space-saving can	SWBA-2-1L-FSC SWB-2-1L-FSC SW55-2-2L-FSC SW60-2-2L-FSC SW90-2-2L-FSC SWBA-2-2L-FSC SWBA-2-2L-FSC	SW55-2-1L-FSC
Funnel with hinged lid with 1 activated charcoal cartridge, 2 inlets 1/8" (3.2mm) and 2 leak fittings (6 to 9mm)	Space-saving can S55 S60 S90 Metallic barrel space-saving can	SWB-2-1L-FSC SW55-2-2L-FSC SW60-2-2L-FSC SW90-2-2L-FSC SWBA-2-2L-FSC SWB-2-2L-FSC	SW55-2-1L-FSC
Funnel with hinged lid with 1 activated charcoal cartridge, 2 inlets 1/8" (3.2mm) and 2 leak fittings (6 to 9mm) Sp Funnel with hinged lid with 1 activated charcoal cartridge,	S55 S60 S90 Metallic barrel	SW55-2-2L-FSC SW60-2-2L-FSC SW90-2-2L-FSC SWBA-2-2L-FSC SWB-2-2L-FSC	
Activated charcoal cartridge, 2 inlets 1/8" (3.2mm) and 2 leak fittings (6 to 9mm) Funnel with hinged lid with 1 activated charcoal cartridge,	S60 S90 Metallic barrel	SW60-2-2L-FSC SW90-2-2L-FSC SWBA-2-2L-FSC SWB-2-2L-FSC	
Activated charcoal cartridge, 2 inlets 1/8" (3.2mm) and 2 leak fittings (6 to 9mm) Funnel with hinged lid with 1 activated charcoal cartridge,	S90 Metallic barrel	SW90-2-2L-FSC SWBA-2-2L-FSC SWB-2-2L-FSC	
2 inlets 1/8" (3.2mm) and 2 leak fittings (6 to 9mm)	Metallic barrel	SWBA-2-2L-FSC SWB-2-2L-FSC	
Funnel with hinged lid with 1 activated charcoal cartridge,	pace-saving can	SWB-2-2L-FSC	
Funnel with hinged lid with 1 activated charcoal cartridge,			
activated charcoal cartridge,	S55		
activated charcoal cartridge,		SW55-2-3L-FSC	
•	S60	SW60-2-3L-FSC	
2 inlets 1/8" (3,2mm) and 3 leak fittings (6 to 9mm)	S90	SW90-2-3L-FSC	
, ,	Metallic barrel	SWBA-2-3L-FSC	
Sp	pace-saving can	SWB-2-3L-FSC	
Create your own reference ! SW55 - 2 - 2L - FSC Type of thread 2 leak fittings (6 to 9 mm) 2 inlets Funnel with hinged lid with 1 activated charcoal cartridge		(exemple: 2	cial configuration entries of different contact us directly!

FUNNEL WITH STOPCOCK


It secures the storage of your residual solvents, equipped with an activated charcoal cartridge for filtration of unobstructed vapors. It is possible to add entries and / or leak fittings by creating your own reference (see below).

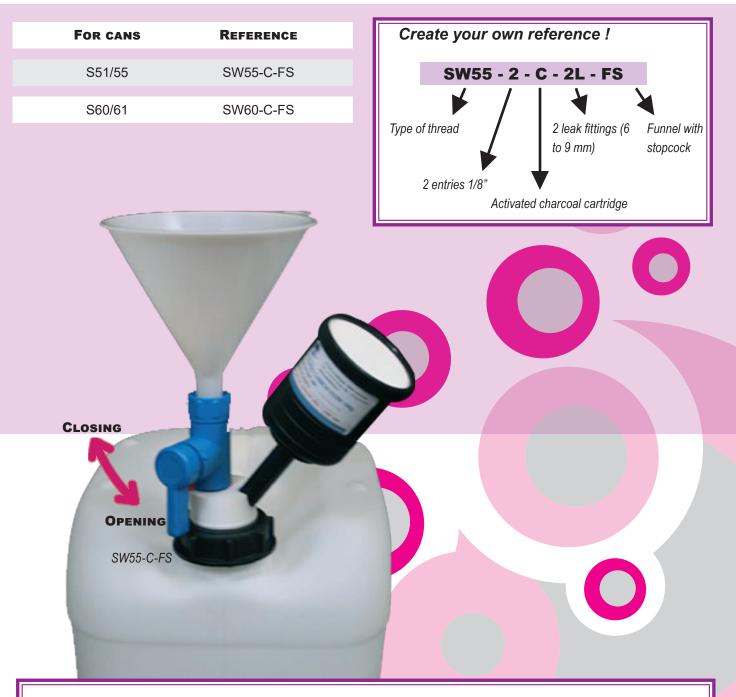

<u>Note</u>: Standard references are not equipped with additional inlets. For another configuration, follow the explanation above and create your own reference!

FUNNEL WITH BALL VALVE WITH HIGH SPILL CAPACITY

Funnel with ball valve allowing to quickly pour used solvents, for plastic can from 5 to 50L.

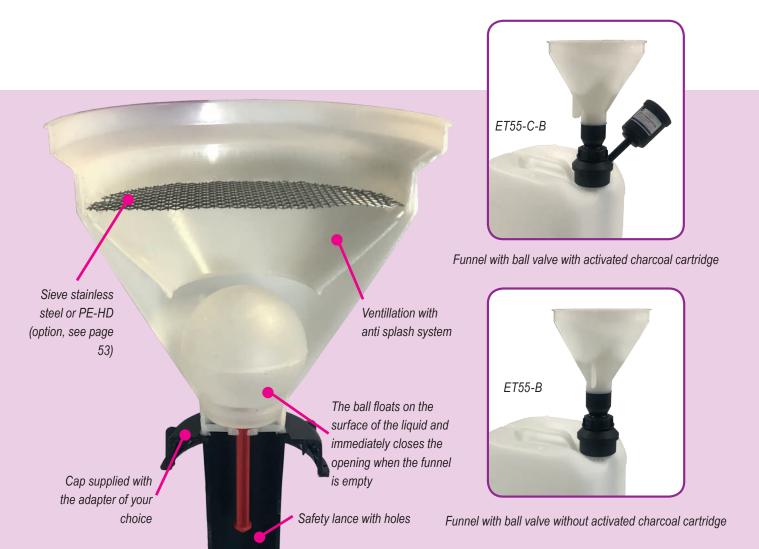


	REFERENCE WITHOUT	REFERENCE WITH
DESCRIPTION	CARTRIDGE	CARTRIDGE *
Funnel with ball valve for can S55	ET55-B	ET55-C-B
Funnel with ball valve for can S60	ET60-B	ET60-C-B
Funnel with ball valve for can S90	ET90-B	ET90-C-B
Funnel with ball valve for space-saving can S50	ETB-B	ETB-C-B
Funnel with ball valve for metallic barrel	ETBA-B	ETBA-C-B

* Activated charcoal cartridge are to be ordered separately http://www.smartcaps.fr

FUNNEL WITH BALL VALVE WITH SMART WASTE CAP FONCTION

Funnel with ball valve for 5 to 20L can and metallic barrel, with or without external solvent inlet. It can be equipped with the same functions as a conventional Smart Waste Cap (leak fittings, activated charcoal cartridge, inlets) to directly receive the effluent from your HPLC chains without passing through the funnel.

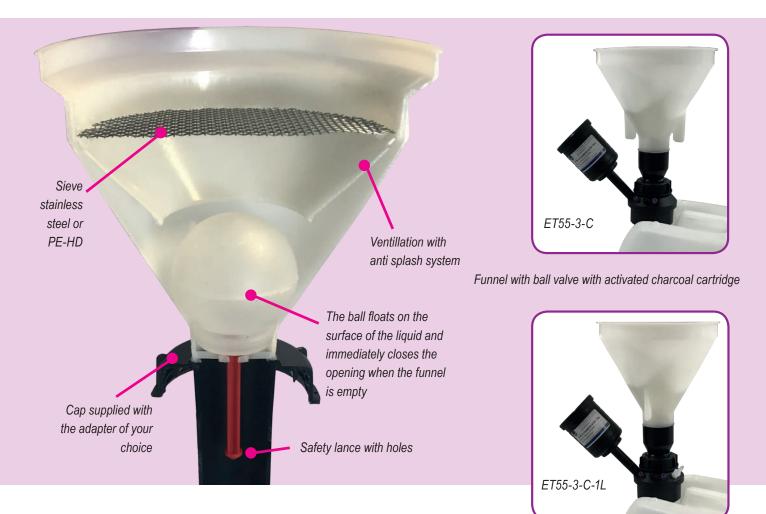

	For any other configuration, contact us		
E	DESCRIPTION	THREAD	Reference
		S55	ET55-2-C-1L
	Funnel with ball valve with an adapter for the activated	S60	ET60-2-C-1L
	charcoal cartridge, 2 fittings 1/8" and 1 leak fitting (6 to 9 mm)	S90	ET90-2-C-1L
		NALGENE® B83	ETB83-2-C-1L
52	A.I.T. FRANCE		http://www.smartcaps.fr

FUNNEL WITH BALL VALVE WITH HIGH SPILL CAPACITY FOR BARREL



DESCRIPTION	Reference
Sieve in polyhethylene (HD-PE)	TAM-PE *
Sieve in stainless steel	TAM-LUX *

* Add the option reference at the end of the funnel reference.





ELECTRICALLY CONDUCTIVE FUNNEL WITH BALL VALVE



DESCRIPTION	THREAD	REFERENCE
	S55	ET55-B-EL
Electrically conductive funnel with ball valve	S60	ET60-B-EL
equipped with a stainless steel sieve	S90	ET90-B-EL
	NALGENE® B83	ETB83-B-EL
	S55	ET55-C-EL
Electrically conductive funnel with ball valve equipped with	S60	ET60-C-EL
a stainless steel sieve and an adapter for an activated charcoal cartridge	S90	ET90-C-EL
	NALGENE [®] B83	ETB83-C-EL


ELECTRO CONDUCTIVE CANS



VOLUME (LITER)	THREAD	REFERENCE
10	S60 / DIN 61	B10-61-EL
20	S60 / DIN 61	B20-61-EL
30	S60 / DIN 61	B30-61-EL

Notes

••••	 •••••••		
• • • • • • •	 	• • • • • • • • • • • • • • • • • • • •	
•••••	 		
• • • • • • •	 • • • • • • • • • • • • • • • •		
•••••	 		
• • • • • • •	 		
• • • • • • •	 		

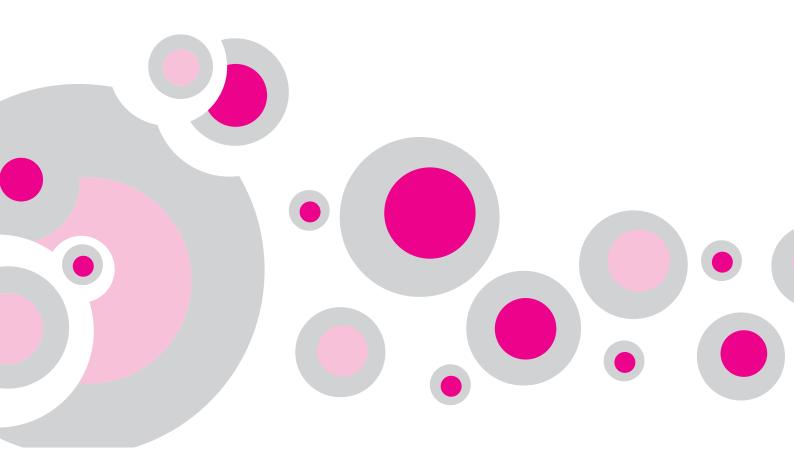

RETENTION TRAYS & CANS CANS IN HD-PE CHIMICALLY INERT

Volume (Liter)	THREAD	REFERENCE
2,5 L	GL45	B2-45
5 L	GL45	B5-45
5 L	S55	B5-55
10 L	S55	B10-55
20 L	S60	B20-60
5 L (space-saving can)	S50	BS-GP





RETENTION TRAYS & CANS


RETENTION **T**RAYS

Avoid accidental overflow. A.I.T offers 3 sizes of retention trays. With a grid to keep the can dry.



	DIMENSION	(мм)			
DESCRIPTION			CAPACITY	Reference	
	EXTERNAL	INTERNAL			
Little Retention Tray	H.70 × D.347 × W.230	D.280 × W.165	3L	BAC-R-3	
Medium Retention Tray	H.75 × D.437 × W.285	D.380 × W.220	5L	BAC-R-5	
Large Retention Tray	H.75 × D.485 × W.335	D.420 × W.265	8L	BAC-R-8	
H : height W : width D : depth					

SPACE-SAVING CAN IN HD-PE

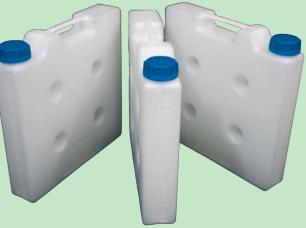

BS-GP

VOLUME (LITER)	THREAD
5	DIN50
DIMENSION (MM) EXT.	Reference

Reference	Page	Reference	Page	Reference	Page
AD-2T-GL45	16	B20-61-EL	54	ET60-B	51
AD-32-45-P	18	B2-45	56	ET60-B-EL	54
AD-32-45-T	18	B30-61-EL	54	ET60-C-B	51
AD-38-45-T	18	B5-45	56	ET60-C-EL	54
AD-38-45-P	18	B5-55	56	ET90-2-C-1L	52
AD-40-32-P	38	BAC-R-3	57	ЕТ90-В	51
AD-40-38-P	38	BAC-R-5	57	ET90-B-EL	54
AD-40-45-P	18	BAC-R-8	57	ET90-C-B	51
AD-40-45-P	38	BS-GP	57	ET90-C-EL	54
AD-40-45-T	18	CC-ONE-1BL	21, 39	ETB83-2-C-1L	52
AD-45-32-P	38	CC-ONE-1V	21, 39	ETB83-B	51
AD-45-38-P	38	CC-ONE-2BL	21, 39	ETB83-B-EL	54
AD-45-55-P	38	CC-ONE-2V	21, 39	ETB83-C	51
AD-45-EM	10	CR-1	42	ETB83-C-EL	54
AD45-PAIL-1L	36	CR-3	42	ETBA-1L	53
AD-4T-GL45	16	CR-6	42	ETBA-2L	53
AD-55-60-P	38	DIS 10	11	ETBA-B	51, 53
AD-55-90-P	38	DIS 25	11	ETBA-C-B	51
AD-60-90-P	38	ET55-2-C-1L	52	ETB-B	51
AD-CPC-JUST	27	ET55-B	51	ETB-C-B	51
AD-D-JUST	27	ET55-B-EL	54	F16-1	22
B10-55	56	ET55-C-B	51	F14	22
B10-61-EL	54	ET55-C-EL	54	F16	22
B20-60	56	ET60-2-C-1L	52	F18	22

Reference	Page	Reference	Page	Reference	Page
F18-25	22	GL5000-C	19	QC32-18-CD-C	45
F316	22	GL500-A	19	QC64-14-CC	45
FC-100	24	GL500-C	19	QC64-14-CC-C	45
FC-100-C	24	KEY-55-60	56	QC64-14-CD	45
FC-100-KOH	24	MULTI-MOL-6-4L	36	QC64-14-CD-C	45
FC-180	26	QC32-116-CD	45	QC64-38-CC	45
FC-25	24	QC32-116-CD-C	45	QC64-38-CC-C	45
FC-250	25	QC32-14-CC	45	QC64-38-CD	45
FC-25-C	24	QC32-14-CC-C	45	QC64-38-CD	45
FC-50	24	QC32-14-CD	45	QC64-38-CD-C	45
FC-500	25	QC32-14-CD-C	45	R14	22
FC-50-C	24	QC32-14-TFP	45	R16	22
FC-90	26	QC32-14-TFP-C	45		
FC-AD-CPC	27	QC32-18-CD	45		
FC-AD-JUST	27				
FC-D	26				
FC-D-90	26				
GL1000-A	19				
GL1000-C	19				
GL2000-A	19				
GL2000-C	19				
GL250-A	19				
GL250-C	19				
GL5000-A	19				
http://www.smartcaps.fr				A.I.T. FRA	NCE 59

Reference	Page	Reference	Page	
R18	22	SHCB53-FSP-B	15	
R316	22	SHCBA-FSP-A	15	
RC-34	41	SHCBA-FSP-C	15	
RC-56	41	SHCBA-FSP-B	15	
RC-69	41	SHC-1-EM	10	
RC-911	41	SHC-P-A	13	
RCY-56	41	SHC-P-B	13	
RCY-69	41	SHC-1PV	11	
RD-34	41	SHC-1-SEPT	12	
RD-56	41	SHC-2	9	
RD-69	41	SHC-2-EM	10	
RD-911	41	SHC-P-A-C	13	
SHC-1	9	SHC-P-2B	13	
SHC45-FSP-A	15	SHC-2PV	11	
SHC45-FSP-C	15	SHC-2-SEPT	12	
SHC45-FSP-B	15	SHC-3	9	
SHC55-FSP-A	15	SHC-3-EM	10	
SHC55-FSP-C	15	SHC-3PV		
SHC55-FSP-B	15	SHC-SV-PREP	20	
SHC60-FSP-A	15			
SHC60-FSP-C	15			
SHC60-FSP-B	15			
SHCB53-FSP-A	15			
SHCB53-FSP-C	15			
60 A.I.T. FRANCE				http://www.smartc

60 A.I.T. FRANCE

Reference	Page	Reference	Page	Reference	Page
SHC-3-SEPT	12	SHC-V4	10	SW55-FSC	48
SHC-4	9	SW45-1-C-1L	30	SW60-1-C-1L	31
SHC-4-EM	10	SW45-1-C-2L	30	SW60-1-C-4L	31
SHC-4PV	11	SW45-2-C-1L-LI	35	SW60-2-1L-FSC	49
SHC-5	9	SW45-2-C	30	SW60-2-2L-FSC	49
SHC55-P-A	14	SW45-2-C-2L	30	SW60-2-3L-FSC	49
SHC55-P-A-B	14	SW45-4-C-1L	30	SW60-2-C	31
SHC55-P-A-C	14	SW45-5-C	30	SW60-2-C-1L	31
SHC55-P-B	14	SW45-C-2L	30	SW60-2-C-4L	31
SHC55-P-2A-B	14	SW45-C-3L	30	SW60-5-C-1L	31
SHC55-P-3A	14	SW55-1-C-3L	29	SW60-6-C	31
SHC60-P-A	14	SW55-2-1L-FSC	49	SW60-C-1L	31
SHC60-P-A-C	14	SW55-2-2L-FSC	49	SW60-C-FS	50
SHC60-P-B	14	SW55-2-3L-FSC	49	SW60-FSB	48
SHC60-P-2A-C	14	SW55-2-C-1L	29	SW60-FSC	48
SHC-SV	20	SW55-2-C-2L	29	SW90-1-C-1L-Z	44
SHC-SV10	20	SW55-3-C-2L	29	SW90-2-1L-FSC	49
SHC-SV50	20	SW55-3-C-3L	29	SW90-2-2L-FSC	49
SHC-SV-AB	20	SW55-4-C	29	SW90-2-3L-FSC	49
SHC-SV-AB10	20	SW55-4-C-1L	29	SW90-2-C	33
SHC-SV-AB50	20	SW55-5-C	29	SW90-2-C-1L	33
SHC-V1	10	SW55-C-Z	44	SW90-2-C-2L	33
SHC-V2	10	SW55-C-FS	50	SW90-3-C	33
SHC-V3	10	SW55-FSB	48	SW90-3-C-1L	33

Reference	Page	Reference	Page	
SW90-3-C-2L	33	SWBA-2-3L-FSC	49	
SW90-4-C	33	SWBA-FSB	48	
SW90-4-C-1L	33	SWBA-FSC	48	
SW90-FSB	48	SWB-FSB	48	
SW90-FSC	48	SWB-FSC	48	
SWB-1-C-1L	31	TAM-LUX	53	
SWB-1-C-3L	31	TAM-PE	53	
SWB-2-1L-FSC	49	TER	54	
SWB-2-2L-FSC	49	TP16	23, 43	(
SWB-2-3L-FSC	49	TP14	23, 43	
SWB-2-C-1L	31	TP18	23, 43	
SWB-3-C	31	TP316	23, 43	
SWB-5-C	31	TYG-127	23, 43	
SWB-5-C-1L	31	TYG-64	23, 43	
SWB53-2-C	34	VP-CA	43	
SWB53-4-C-1L	34	VP-F18	22, 43	
SWB53-C-1L	34	VP-LEAK	43	
SWB53-C-3L	34	Y-11	42	
SWB83-2-C	34	Y-3	42	
SWB83-4-C-1L	34	Y-5	42	
SWB83-C-1L	34	Y-7	42	
SWB83-C-3L	34	Y-9	42	
SWBA-2-1L-FSC	49			
SWBA-2-2L-FSC	49			



TO PROTECT YOUR HEALTH



DISTRIBUTED BY:

Greyhound Chromatography & Allied Chemicals 6 Kelvin Park, Birkenhead, Merseyside, CH41 1LT, UK Tel: (+44) 0 151 649 4000 • Fax: (+44) 0 151 649 4001 Email: info@greyhoundchrom.com • Web: www.greyhoundchrom.com